

PURITAN BENNETT

840™ VENTILATOR SYSTEM



The Puritan Bennett® 840™ ventilator system is the flagship product in our line of critical care ventilators. It is highly responsive and offers superior comfort, delivering sensitive, precise breaths to critically ill neonatal through adult patients.

Be Certain.

WHY THE 840™ VENTILATOR IS THE WORLD'S BEST-SELLING CRITICAL CARE VENTILATOR



Patient comfort

With highly responsive proportional solenoid valves, an active exhalation valve and state-of-the-art flow sensors, the 840 ventilator provides exceptional responsiveness and control. These features deliver sensitive and precise breaths, and in turn, superior comfort to patients.

Patient and clinician safety

SmartAlert™ Alarm System: Prioritized alarm annunciation distinguishes primary alarms from secondary, “dependent” alarms.

- System indicates the urgency level of each alarm to help you efficiently respond to alarm conditions.

Protection from contaminants: A special filter heating process shields patients and clinicians from exposure to viruses and bacteria from exhaled gases.

- Expiratory filter is capable of filtering viruses such as SARS.



Circuit disconnect alarm: Special circuit disconnect alarm does not rely on low pressure or tidal volume to activate.

- When a circuit disconnects, gas flow stops to prevent the spread of contaminants.
- Technicians are alerted so they can quickly identify the source of the disconnection.

Self-diagnostic testing: Subsystems are thoroughly checked 100 times per second.

- Message display notifies technicians of potential failures if any drift in performance occurs while ventilating.

Seamless communication

The 840 ventilator supports communication with all major patient monitoring and hospital information systems. It also provides seamless recorded data transfer into a patient's electronic medical record.



Puritan Bennett Clinivision Mobile Patient Charting System: The 840 ventilator teams up with the Clinivision system to automate respiratory care patient charting and help balance RT workloads. The Clinivision system:

- Offers the latest in wireless, mobile clinical information technology.
- Puts patient information right where you need it—at your fingertips.
- Produces more accurate records and improved charge capture.
- Supports evidence-based medicine for improved patient care.


Ease of use

The innovative design of the 840 ventilator makes it easy for technicians to set up and navigate. It minimizes training requirements for your staff.

- DualView™ LCD touchscreens display monitored data separately from ventilator settings, so clinicians can change settings and view the impact in one easy place.
- Short self-test feature resets the 840 ventilator settings in between patient use in about two minutes. The test checks for proper ventilation operation and calculates circuit compliance.



You can change display selections and scaling even while waveforms are frozen.

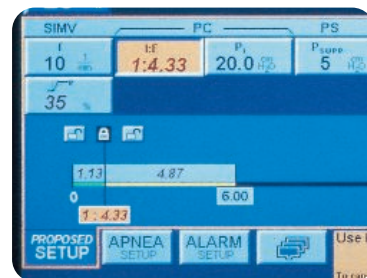

Pressure-based breath delivery allows you to select inspiratory time, I:E ratio, or expiratory time as a direct setting.



Shown here is the PAV+ software option, a spontaneous breath type that adjusts ventilator support based on patient demand.

For more information on the 840 ventilator, please contact your Puritan Bennett sales representative or call Customer Service at 1-800-635-5267.

www.puritanbennett.com

840 Ventilator System

What our customers are saying about the 840 ventilator...

"After much comparison and study, with input from our staff, nurses and physicians, we chose the Puritan Bennett 840 ventilator to replace our aging fleet of 7200s. It offered the most technology coupled with ease of use and affordable price. In addition, the education, follow-up and service has been superb. We purchased 25 ventilators and are in the process of placing an order for 25 more."

– JAN R. GODDARD, MSM, RRT

"The Puritan Bennett 840 ventilator is user-friendly, clinically versatile, and gives us the ability to offer our patients excellent care. We depend on Puritan Bennett's proven record for dependable "work horse" ventilators as well as for their superior service."

– SUE HENRY, RRT



tyco
Healthcare

**Puritan
Bennett**

4280 Hacienda Drive
Pleasanton, CA 94588
Tel 925.463.4000
Toll Free 1.800.635.5267
www.puritanbennett.com

Tyco Healthcare UK LTD.
154 Fareham Road
Gosport, UK PO13 0AS
Tel +44.1329.224000

Specifications subject to
change without notice.

© 2006 Nellcor Puritan
Bennett Incorporated.
All rights reserved.

*Proportional Assist and PAV
are registered trademarks of
The University of Manitoba
and are used under license
by Puritan Bennett. All other
trademarks belong to Tyco
Healthcare Group LP or an
affiliate.

O.c 2070v3-0506
VE03200